**SIEMENS**

# Orion™ 2 RIC10

## User Guide

[usa.bestsoundtechnology.com](http://usa.bestsoundtechnology.com)

**Life sounds brilliant.**

## Welcome

### Enhance your Siemens experience.

Learn about your hearing aids and accessories you may find helpful, access the latest information about hearing health, and receive promotional offers by joining the Siemens hearing health community at\* [usa.bestsoundtechnology.com/join](https://usa.bestsoundtechnology.com/join) or:

Learn more  
about Siemens  
products by  
scanning the  
QR code

Thank you for choosing our hearing instruments to accompany you through everyday life. Like anything new, it may take you a little time to become familiar with them.

This guide, along with support from your Hearing Care Professional, will help you understand the advantages and greater quality of life your hearing instruments offer.

To receive the most benefit from your hearing instruments, it is recommended that you wear them every day, all day. This will help you get adjusted to them.

#### CAUTION

It is important to read this user guide and the safety manual thoroughly and completely. Follow the safety information to avoid damage or injury.

\*We respect your privacy. We will not sell, rent, or give your name or e-mail address to anyone. At any time, you can select a link at the bottom of every email to unsubscribe or update profile.

# Content

Welcome	1
Your hearing instruments	4
Instrument type	4
Getting to know your hearing instruments	4
Components and names	5
Controls	7
Settings	8
Batteries	9
Battery size and handling tips	9
Replacing batteries	10
Daily use	12
Turning on and off	12
Inserting and removing the hearing instruments	13
Adjusting the volume	16
Changing the hearing program	16
Special listening situations	17
On the phone	17
Maintenance and care	18
Hearing instruments	18
Ear pieces	19
Professional maintenance	21

Further information	22
Safety information	22
Accessories	22
Troubleshooting	22
Service information	24

## Symbols used in this document

Points out a situation that could lead to serious, moderate, or minor injuries.

Indicates possible property damage.

Advice and tips on how to handle your device better.

# Your hearing instruments

**i** This user guide describes optional features that your hearing instruments may or may not have. Ask your Hearing Care Professional to indicate the features that are valid for your hearing instruments.

## Instrument type

Your hearing instruments are a RIC (Receiver-in-Canal) model. The receiver is placed within the ear canal and connected to the instrument via a receiver cable. The instruments are not intended for infants, small children or mentally disabled persons.

## Getting to know your hearing instruments

It is recommended to familiarize yourself with your new hearing instruments. With the instruments in your hand try using the controls and note their location on the instrument. This will make it easier for you to feel and press the controls while wearing the hearing instruments.

## Components and names

- |                          |   |
|--------------------------|---|
| 1 Ear piece              | 6 Battery compartment                                     |
| 2 Receiver               | 7 Side indicator<br>(red = right ear,<br>blue = left ear) |
| 3 Receiver cable         | 8 Receiver connection                                     |
| 4 Microphone openings    |   |
| 5 Push button (optional) |   |

You can use standard or custom-made ear pieces.

Standard ear pieces		Size
<input type="checkbox"/>		Click Dome single (open or closed)
<input type="checkbox"/>		Click Dome semi-open
<input type="checkbox"/>		Click Dome double


You can easily exchange the standard ear pieces. Read more in section "Maintenance and care".

Custom-made ear pieces	
<input type="checkbox"/>	 Custom shell
<input type="checkbox"/>	 Click Mold

## Controls

With the controls you can, for example, switch hearing programs. Your hearing instruments either have a push button or no control.

Your Hearing Care Professional has programmed your desired functions to the push button.

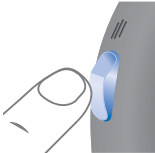Control	Left	Right
Push button	<input type="checkbox"/>	<input type="checkbox"/>
		
No control	<input type="checkbox"/>	<input type="checkbox"/>

**i** You can use a smartphone app for remote control.

Push button function	L	R
Press <1 second:		
Program change	<input type="checkbox"/>	<input type="checkbox"/>
Press >3 seconds:		
Standby/turn on	<input type="checkbox"/>	<input type="checkbox"/>

L = Left, R = Right

## Settings

Hearing programs
1
2
3
4
Read more in section "Changing the hearing program".
Features
<input type="checkbox"/> <b>Power-on delay</b> enables whistle-free insertion of hearing instruments. Read more in section "Turning on and off".
<input type="checkbox"/> The <b>tinnitus therapy feature</b> generates a soft noise to distract you from your tinnitus.

## Batteries

When the battery is low the sound becomes weaker or you will hear an alert signal. The battery type will determine how long you have to replace the battery.

### Battery size and handling tips

Ask your Hearing Care Professional for recommended batteries.

Battery size:	10
<ul style="list-style-type: none"><li>■ Always use the correct battery size for your hearing instruments.</li><li>■ Remove the batteries if you intend not to use the hearing instruments for several days.</li><li>■ Always carry spare batteries.</li><li>■ Remove empty batteries immediately and follow your local recycling rules for battery disposal.</li></ul>	


## Replacing batteries

Removing the battery:

- ▶ Open the battery compartment.

- ▶ Use the magnet stick to pull out the battery. The magnet stick is available as an accessory.

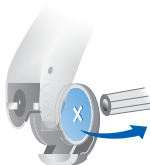Inserting the battery:

- ▶ If the battery has a protective film, remove it only when you are ready to use the battery.

- ▶ Insert the battery with the "+" symbol facing up (refer to the picture).

- ▶ Carefully close the battery compartment. If you feel resistance, the battery is not inserted correctly.

Do not attempt to close the battery compartment by force. It can be damaged.

# Daily use

## Turning on and off

You have the following options to turn your hearing instruments on or off.

Via battery compartment:

- ▶ Turning on: Close the battery compartment.  
The **default** volume and hearing program are set.
- ▶ Turning off: Open the battery compartment to the first catch.

Via push button:

- ▶ Turning on or off: Press and hold the push button for several seconds. Refer to section "Controls" for the settings of your controls.

After turning on, the **previously** used volume and hearing program remain set.

When wearing the hearing instruments, an alert tone can indicate when an instrument is being turned on or off.

When the **power-on delay** is activated, the hearing instruments turn on after a delay of several seconds. During this time you can insert the hearing instruments into your ears without experiencing unpleasant feedback whistling.

The "power-on delay" can be activated by your Hearing Care Professional.

## Inserting and removing the hearing instruments

Your hearing instruments have been fine-tuned for your right and left ear. Colored markers indicate the side:

- red marker = right ear
- blue marker = left ear


Inserting a hearing instrument:

- ▶ Hold the receiver cable at the bend closer to the ear piece.
- ▶ Carefully push the ear piece in the ear canal **1**.
- ▶ Twist it slightly until it sits well.

Open and close your mouth to avoid accumulation of air in the ear canal.

- ▶ Lift the hearing instrument and slide it over the top of your ear **2**.
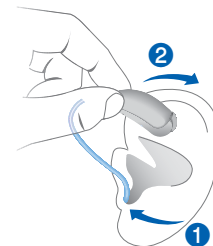
### CAUTION

Risk of injury!


- ▶ Insert the ear piece carefully and not too deeply into the ear.

# i

- It may be helpful to insert the **right** hearing instrument with the **right** hand and the **left** hearing instrument with the **left** hand.
- If you have problems inserting the ear piece, use the other hand to gently pull your earlobe downwards. This opens the ear canal and eases insertion of the ear piece.


The optional retention strand helps to securely retain the ear piece in your ear. To position the retention strand:

- ▶ Bend the retention strand and position it carefully into the bottom of the bowl of your ear (refer to the picture).

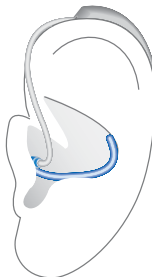Removing a hearing instrument:

- ▶ Lift the hearing instrument and slide it over the top of your ear ①.
- ▶ If your hearing instrument is equipped with a custom shell or a Click Mold, remove it by pulling the small removal cord toward the back of your head.
- ▶ For all other ear pieces: Grip the receiver in the ear canal with two fingers and pull it carefully out ②.
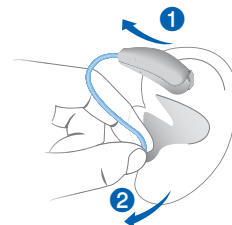
## CAUTION

Risk of injury!

- ▶ In very rare cases the ear piece could remain in your ear when removing the hearing instrument. If this happens, have the ear piece removed by a medical professional.

Clean and dry your hearing instruments after usage. Read more in section "Maintenance and care".

## Special listening situations

### Adjusting the volume

Your hearing instruments automatically adjust the volume to the listening situation.

You can use a smartphone app for volume control.

### Changing the hearing program

Depending on the listening situation, your hearing instruments automatically adjust their sound.


Your hearing instruments may also have several hearing programs which allow you to change the sound, if needed. An optional signal tone can indicate the program change.

- To change the hearing program, press the push button briefly.

Refer to section "Controls" for the setting of your controls. Refer to section "Settings" for a list of your hearing programs.

### On the phone

When you are on the phone, hold the telephone receiver slightly above your ear. The hearing instrument and telephone receiver have to be aligned. Turn the receiver slightly so that the ear is not completely covered.
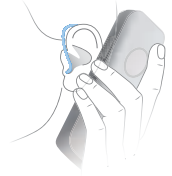
### Telephone program

You may prefer a certain volume when using the phone. Ask your Hearing Care Professional to configure a telephone program.

- Switch to the telephone program whenever you are on the phone.

If a telephone program is configured for your hearing instruments, it is listed in section "Settings".

# Maintenance and care

To prevent damage it is important that you take care of your hearing instruments and follow a few basic rules, which will soon become a part of your daily routine.

## Hearing instruments

### Drying and storage

- ▶ Dry your hearing instruments overnight.
- ▶ Ask your Hearing Care Professional for recommended drying products.
- ▶ For longer periods of non-use, store your hearing instruments with open battery compartment and batteries removed in a drying system to avoid the adverse effects of moisture.

### Cleaning

Your hearing instruments have a protective coating. However, if not cleaned regularly you may damage the hearing instruments or cause personal injury.

- ▶ Clean your hearing instruments daily with a soft, dry tissue.
- ▶ Never use running water or immerse the devices in water.
- ▶ Never apply pressure while cleaning.
- ▶ Ask your Hearing Care Professional for recommended cleaning products, special care sets, or more information on how to keep your hearing instruments in good condition.

## Ear pieces


### Cleaning

Cerumen (ear wax) may accumulate on the ear pieces. This may affect sound quality. Clean the ear pieces daily.

For all types of ear pieces:


- ▶ Clean the ear piece right after removal with a soft and dry tissue.

This prevents that cerumen becomes dry and hard.

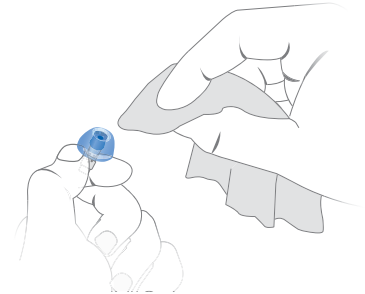For Click Domes:

- ▶ Squeeze the tip of the dome.

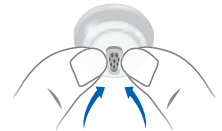
### Exchanging

Exchange standard ear pieces approximately every three months. Replace them sooner if you notice cracks or other changes. The procedure for exchanging standard ear pieces depends on the type of ear piece. In section "Components and names", your Hearing Care Professional has marked your type of ear piece.

### Exchanging Click Domes

If you have problems removing the old Click Dome, use the removal tool or follow the instructions on the packaging of the Click Domes:

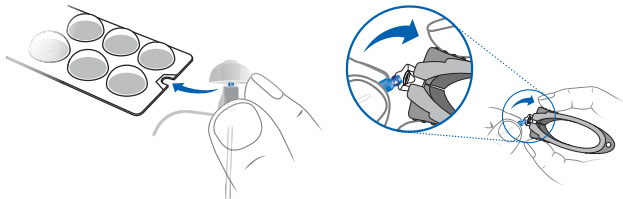The removal tool is available as an accessory.

- Take particular care that the new Click Dome clicks noticeable in place.

### Professional maintenance

Your Hearing Care Professional can perform a thorough professional cleaning and maintenance.

Custom-made ear pieces and wax filters should be exchanged as needed by a professional.

Ask your Hearing Care Professional for your individual recommendation of maintenance intervals and for support.

## Further information

### Safety information

For more safety information, refer to the safety manual supplied with the device.

### Accessories

You can use a smartphone app for remote control.

Consult your Hearing Care Professional for further information.

### Troubleshooting

#### Problem and possible solutions

Sound is weak.

- Replace the empty battery.
- Clean or exchange the ear piece.

Hearing instrument emits whistling sounds.

- Reinsert the ear piece until it fits securely.
- Clean or exchange the ear piece.

Sound is distorted.

- Replace the empty battery.
- Clean or exchange the ear piece.

Hearing instrument emits signal tones.

- Replace the empty battery.

#### Problem and possible solutions

Hearing instrument does not work.

- Turn the hearing instrument on.
- Gently close the battery compartment completely.
- Replace the empty battery.
- Make sure the battery is inserted correctly.
- Power-on delay is active. Wait several seconds and check again.

Consult your Hearing Care Professional if you encounter further problems.

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Information in this brochure is subject to change without notice. Copyright © 2015 Siemens AG. All rights reserved.
La información que contiene este folleto está sujeta a cambios sin aviso. Copyright © 2015 Siemens AG. Todos los derechos reservados.
under Trademark License of Siemens AG
Sivantos, Inc.